

RPA071Hu01 10µg

Recombinant Interleukin 1 Alpha (IL1a)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[**PROPERTIES**]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Ser113~Ala271

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 85%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 350µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.5

Predicted Molecular Mass: 21.2kDa

Accurate Molecular Mass: 21kDa as determined by SDS-PAGE reducing conditions.

[**USAGE**]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[**STORAGE AND STABILITY**]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[**SEQUENCE**]

SAPFSFLS NVKYNFMRII KYEFILNDAL NQSIIRANDQ
 YLTAAALHNL DEAVKFDMDGA YKSSKDDAKI TVILRISKYQ LYVTAQDEDQ
 PVLLKEMPEI PKTITGSETN LLFFWETHGT KNYFTSVAHP NLFIAATKQDY
 WVCLAGGPPS ITDFQILENQ A

[IDENTIFICATION]

ATGCGACCTTTAGCTTCTGAGGATGTGAATCAACTTTATGAGGTCATCAATACGAAATTCCTGAAATGAGGCTCAATCAAAATATAATTCGAGCAATGATCACTACCTCAGCGGCTCCTGCAATGATAATCTGGTGAAGCAATGAATTTGACATGGTCTTTATPAAGTCAATCAAGGATGATGCTAAAATTACCGTG
 S A P F S F L S N V K Y N F M R I I K Y E F I L N D A L N Q S I I R A N D Q Y L T A A A L H N L D E A V K F D M G A Y K S S K D D A K I T V


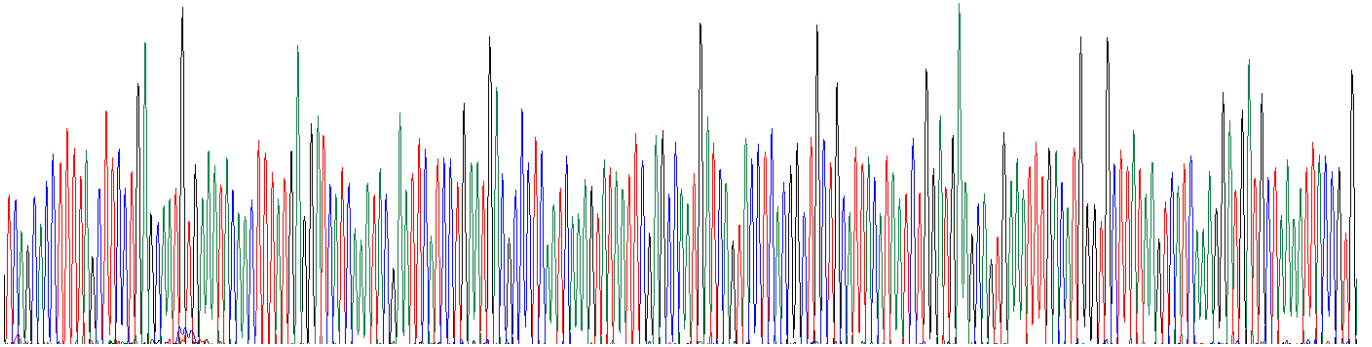


Figure. Gene Sequencing (Extract)

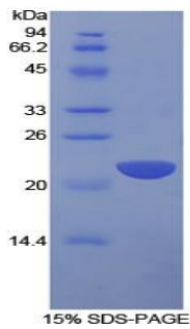

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.